


ViSee® LVM-480 Portable Digital Reading Magnifier

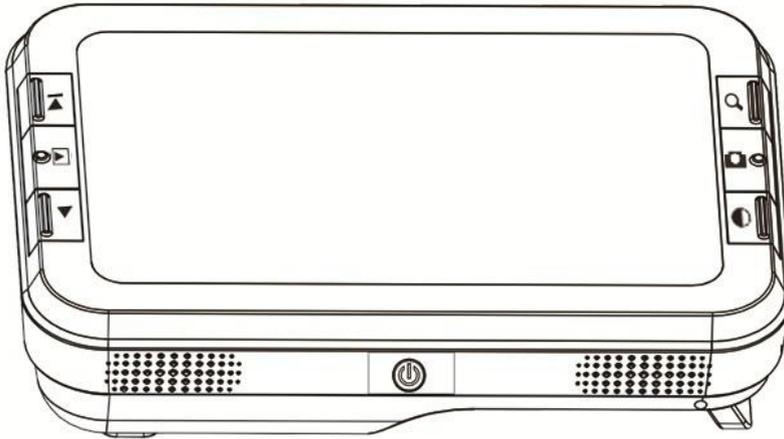

User's Guide

Table of Contents

Precaution	2
Overview	2
Package contents.....	2
Part Names & Functions.....	3
Safety and Maintenance.....	4
Getting Started	4
Loading Battery.....	4
Reading with the Magnifier.....	5
Using the Magnifier	6
Color Mode Selection.....	6
Magnification Adjustment.....	6
Freeze image.....	7
Capture Photo.....	7
Playback.....	7
Delete.....	7
Download to a PC.....	8
Output to a Monitor with AV-in.....	8
Low Battery Indication.....	8
Power-off.....	8
Trouble Shooting	9

Welcome

Thank you for purchasing the ViSee® LVM-480 Portable Digital Magnifier.
Please take a moment to read through this manual. This will help you get the best use of the device.



Precaution

Federal Communications Commission Notice


This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

Overview

Package contains:

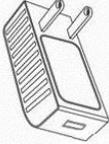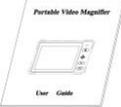		
Portable Digital Magnifier	AC adaptor	Li-ion Battery
		
User's manual	Carrying bag	Cleaning fabric
		
TV-out cable	USB cable	

Part Names & Functions

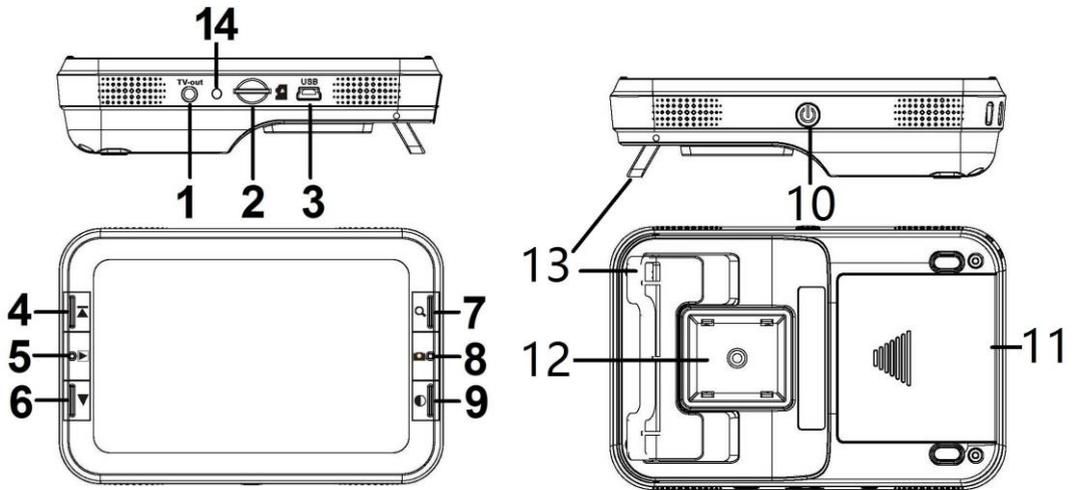

1. TV out
2. Micro-SD slot
3. USB port/DC in
4. Forward/Freeze
5. Playback
6. Backward
7. Magnify

8. Capture/Delete
9. Color mode
10. Power on/off
11. Battery compartment
12. Lens window
13. Foldable stand
14. indicator light

Safety and Maintenance

To keep your Magnifier in good condition, observe the following guidelines:


1. Read the instructions in this user's manual before attempting to use the unit.

2. Do not attempt to repair the Magnifier or open its case. This product contains no user serviceable parts.
3. Keep the Magnifier away from wet areas, rain, and other sources of liquids. Do not submerge the unit.
4. Always operate the Magnifier in a place where the temperature is between 32°F and 113°F.
5. Always store the Magnifier in a cool and dry place at normal room temperature.
6. Do not attempt to force any connector into the unit's USB and Micro-SD ports. All appropriate connectors should fit easily into the ports.

Getting Started

Loading Battery

Please use the lithium battery. Slightly push and pull the battery door latch to open the battery compartment. Load the battery according to the polarities shown inside the battery compartment, and then snap close the battery door.

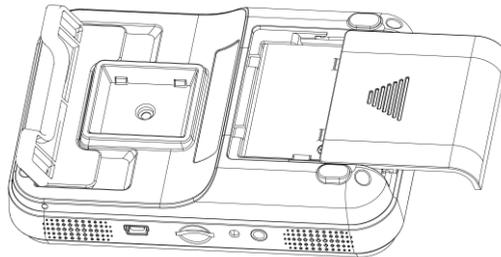


Note: Remove the battery when it is not in use for a long period of time.

Charging the Battery

Before using the Magnifier, you must charge its battery.

1. Connect the adaptor with USB cable. Plug the mini USB port of the power adaptor into the USB port  located on the right side of the unit. Then plug the other end of the power adaptor into a standard AC outlet.

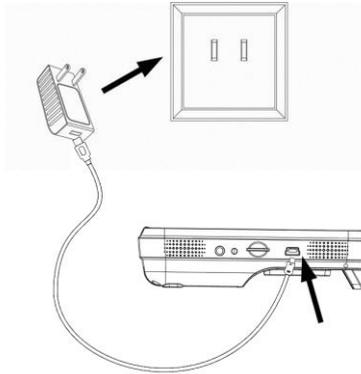

2. When the battery is finished charging, the battery indicator LED will be off. Unplug the unit. Then the battery indicator LED will turn off.
3. You can now use your Magnifier for up to two hours of continuous use. When the battery is low, the low-power indicator  will blink on the right down corner of the screen. Then you have to connect the power adaptor and plug it in to recharge the battery.

Note: Always turn the Magnifier off during charging. The battery will not charge if the Magnifier is left on.

Reading with the Magnifier

To read with the Magnifier, operate as the following:

1. Turn on the unit by pressing the Power on/off switch  located on side.
2. Pull out the foot of the Magnifier until it sits at its maximum range.
3. Place the Magnifier flat against the text you want to read.
4. Slide the unit gradually to the right to read to the end of the line. Slide the unit down and back to the left to read the next line of text.
5. To change the level of magnification, press the Magnify button .
6. To change the colors used on the screen, press the Color Mode button  until you find the mode best suited for your vision needs.

Using the Magnifier

Color Mode selection

Press Color Mode button  to cycle through the Magnifier's 3 different color modes as follows:

- ◆ Full color
- ◆ Enhanced positive image (black on white)
- ◆ Enhanced negative image (white on black)

Magnification Adjustment

Press Magnify button  to cycle through the four level of magnifications
6X, 8X, 10X or 14X

Freeze Image

Press Freeze/Unfreeze  button to switch between freezing and unfreezing the current image.

Capture Photo

Load a micro SD card (not included. the micro SD card is available at many local drug, electronics and department stores and online stores) into the micro SD card slot.

Press snapshot capture  button to capture a photo. The captured photos will be automatically stored in the micro SD card.

Playback

Press Playback button  to view captured photos. Press backward button  and Forward button  to move forward and backward.

Press Playback  button again will return the magnifier to reading mode.

Delete

When press Playback  button to view captured photos and in case you want to delete the pictures, press capture/delete  button to delete it.

Download Captured Pictures to a PC


Connect the Magnifier to a PC with USB Cable to USB port, turn the power on and then you can download the captured pictures to computer.

Output to Monitor with TV-in

Connect the Magnifier with monitor with TV-in by TV video cable, turn the power on and then you can see live video on the monitor.

Low Battery Indication

If battery power is low, there will be a  symbol (blank) on the lower right corner of the screen, indicating that you have to charge the battery.

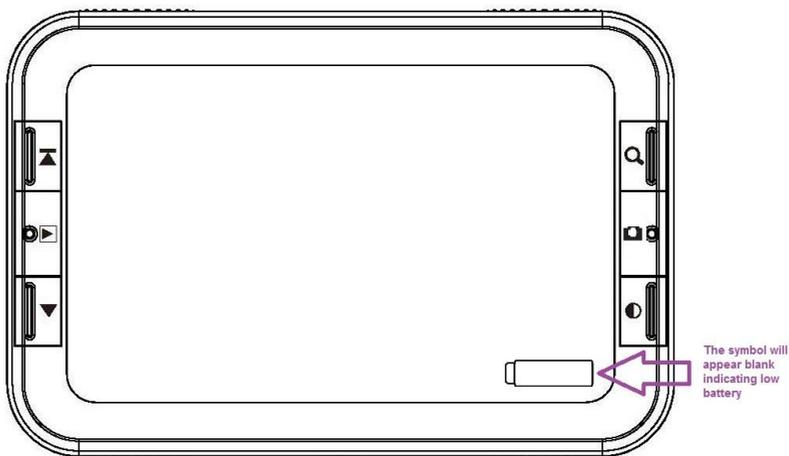

Power Off

Turn off the power of the unit by pressing the Power on/off switch (10) located at the side.

To prolong the battery life, please always turn off the Magnifier when not using it. The Magnifier will turn off automatically if it is not in use for 5 minutes.

Trouble Shooting

Problem	Solution
The screen is blank	Ensure that the power is on. Ensure the battery has enough power. Verify that the item you want to magnify is positioned correctly under the center of the screen
The unit can't be turned on	Change the battery with a new one. Check the installed battery's position
The image on the screen is dark and/or distorted.	Ensure that the unit is positioned flat on the page you are reading. Ensure that the foot of the magnifier flips out to its maximum position.
There are smudges or blemishes on the screen	Clean the screen and the Lens (located on the bottom of the unit).
Text is too small to read	Press the magnify button.

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU.

For more information, please contact the local authority or your retailer where you purchased the product.



ViSee® Magnifiers

ViSee Tech

www.viseetech.com